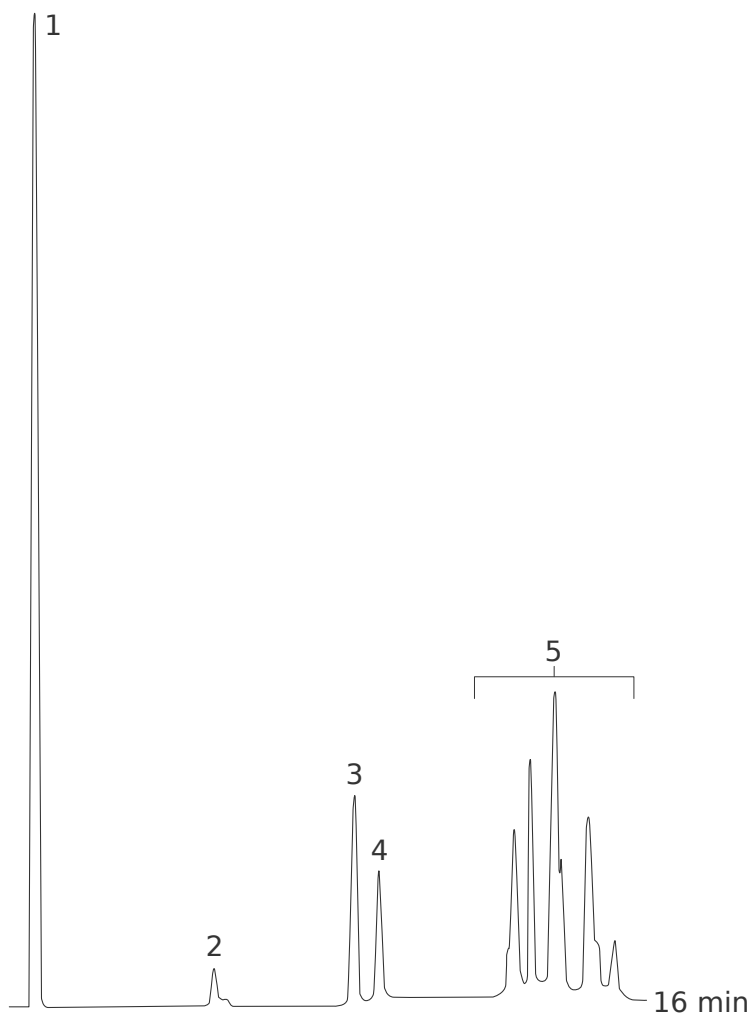




Alkyl Naphthalens

Column: **MEGA-5 - 0.53mm, 1.00µm, 10m**

Catalog Code: S-5-053-100-10

Conditions:

Injection: Split 250°C, 1.0µL.

Detector: FID 250°C.

Oven Temp.: 90°C, 20°C/min, 150°C.

Carrier Gas: Helium, 4.5 mL/min.

Peaks:

1. m-Xylene
2. Naphthalene
3. 2-methylnaphthalene
4. 1-methylnaphthalene
5. others alkyl naphthalens

