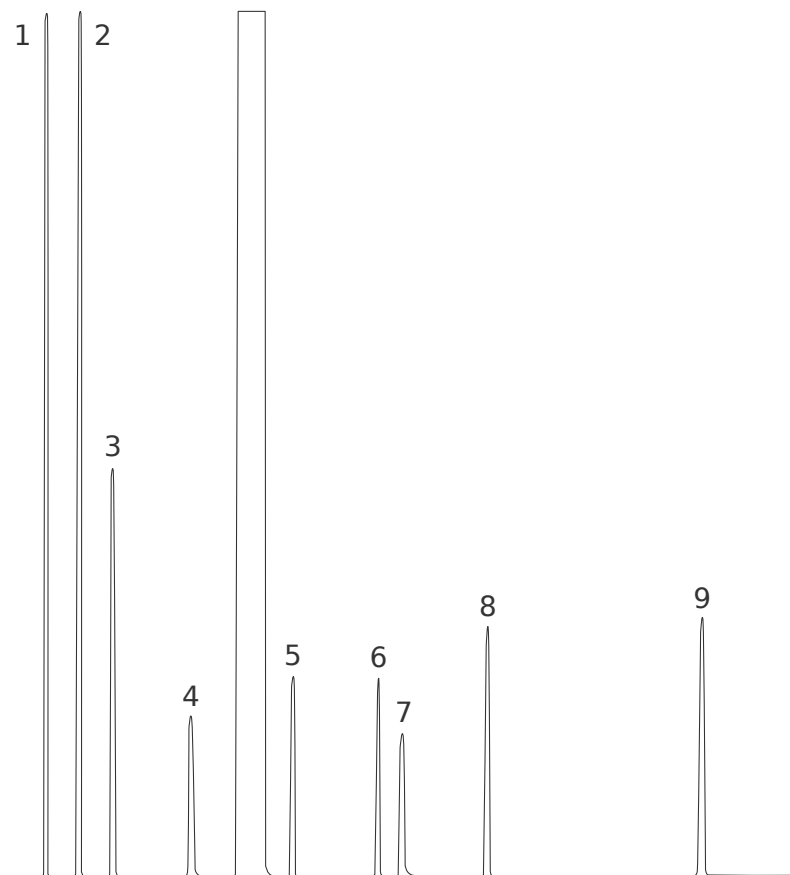




Amines

Column: **MEGA-BASIC - 0.32mm, 0.25 μ m, 25m**

Catalog Code: C-BASIC-032-025-25

Conditions:

Injection: Split 250°C, 1.0 μ L, 1:20 Split Ratio.

Detector: FID 250°C.

Oven Temp.: 35°C, 2.5°C/min, 130°C.

Carrier Gas: Hydrogen, 40kPa.

Peaks:

1. N-ethyldiisopropyl amine
2. Di-N-propyl amine
3. Diallyl amine
4. 2-methyl pyrroline
5. Pyridine
6. 2-methyl diazol
7. Morpholine
8. 4,5-dimethyl diazol
9. Pyrrole

