


Aroclor 1254 - 1260 PCBs (Quick Screening)

Column: **MEGA-SE54 - 0.32mm, 0.25µm, 25m**

Catalog Code: S-SE54-032-025-25

Aroclor 1254

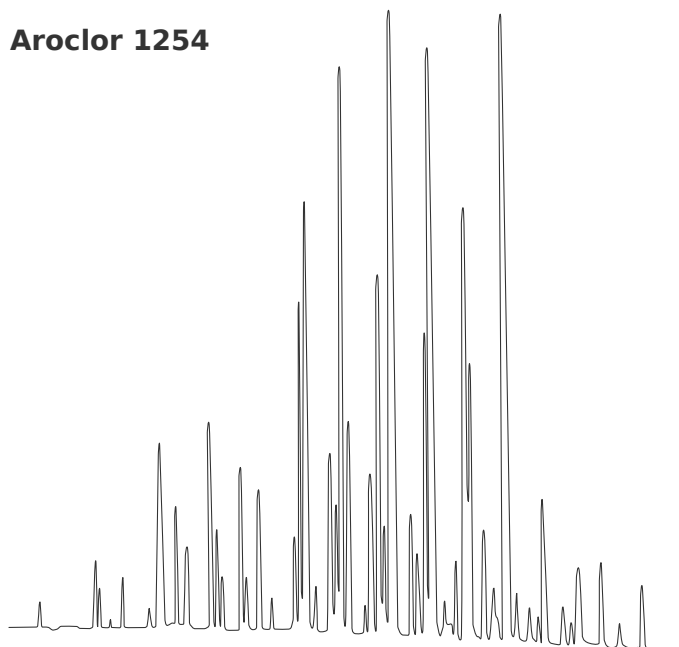

Conditions:

Injection: Split 250°C, 1µL, 1:30 Split Ratio.

Sample: 1000µg/mL in Hexane.

Detector: ECD 300°C, Make up Gas Argon/Methane.

Oven Temp.: 100°C, 3.5°C/min, 250°C.

Carrier Gas: Hydrogen, 70 kPa.

Aroclor 1260

