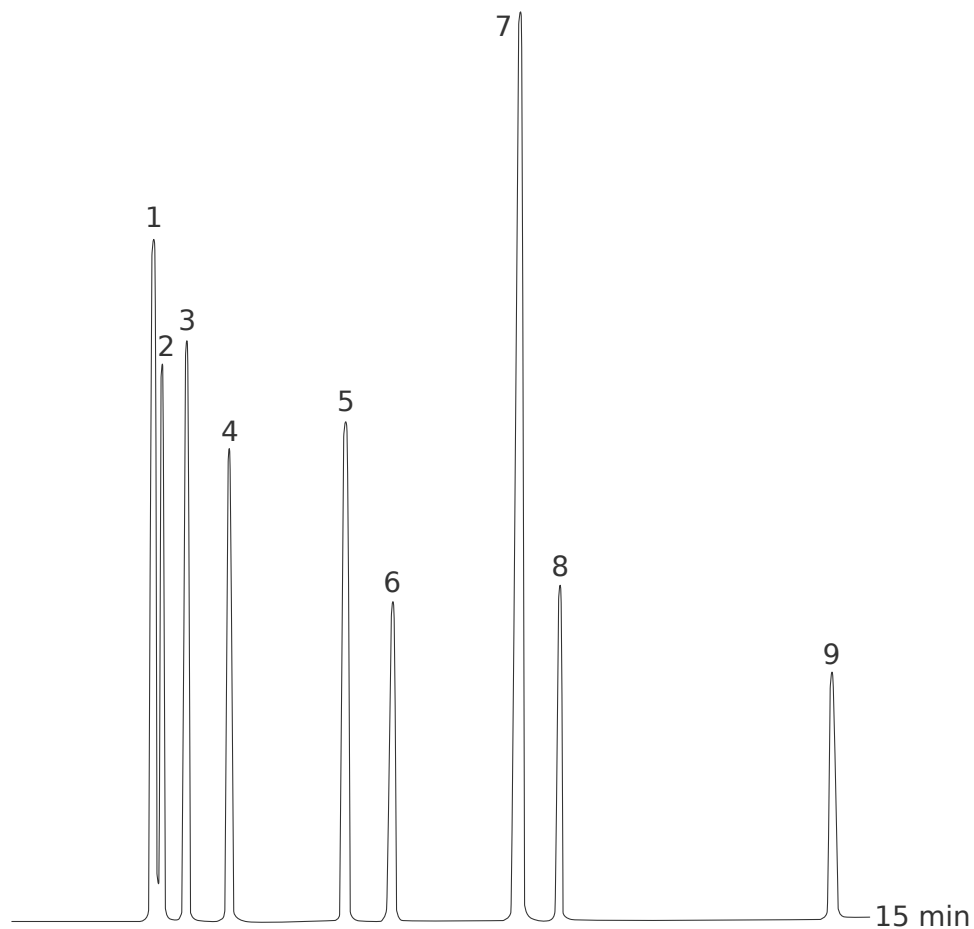




Chlorinated Hydrocarbons EPA Method 612

Column: **MEGA-5 - 0.53mm, 1.50µm, 15m**

Catalog Code: S-5-053-150-15

Conditions:

Injection: On-Column, 1µL.

Sample: 20µg/mL each in Isooctane.

Detector: ECD 250°C.

Oven Temp.: 50°C, 8°C/min, 150°C.

Carrier Gas: Helium, 20 kPa.

Peaks:

1. 1,3-dichlorobenzene
2. 1,4-dichlorobenzene
3. 1,2-dichlorobenzene
4. Hexachloroethane
5. 1,2,4-trichlorobenzene
6. Hexachlorobutadiene
7. Hexachlorocyclopentadiene
8. 2-chloronaphthalene
9. Hexachlorobenzene

