

Delta-Lactones

Chiral Separation on MEGA-DEX Columns

Column: **MEGA-DEX DAC Beta - 0.25mm, 0.25 μ m, 25m**

Catalog Code: DEX-DAC-B-025-025-25

Conditions:


Injection: Split 250°C, 0.5 μ L, 1:50 Split Ratio.

Sample dilution: 1mg/1mL each in n-Hexane.

Contact us for additional information.

Peaks:

1. (+/-) Delta-Octalactone
2. (+/-) Delta-Nonalactone
3. (+/-) Delta-Decalactone
4. (+/-) Delta-Undecalactone
5. (+/-) Delta-Dodecalactone

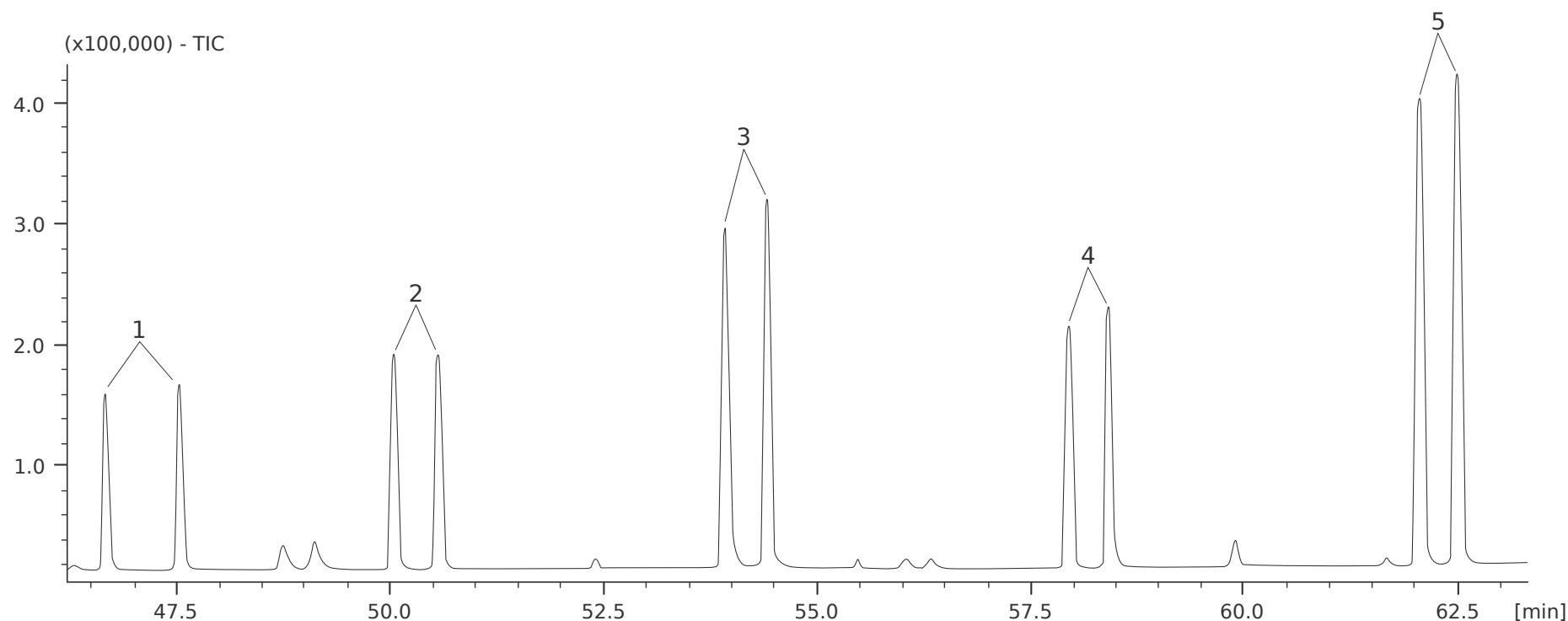

Acknowledgement: Prof. C.Bicchi, C.Cordero et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.

