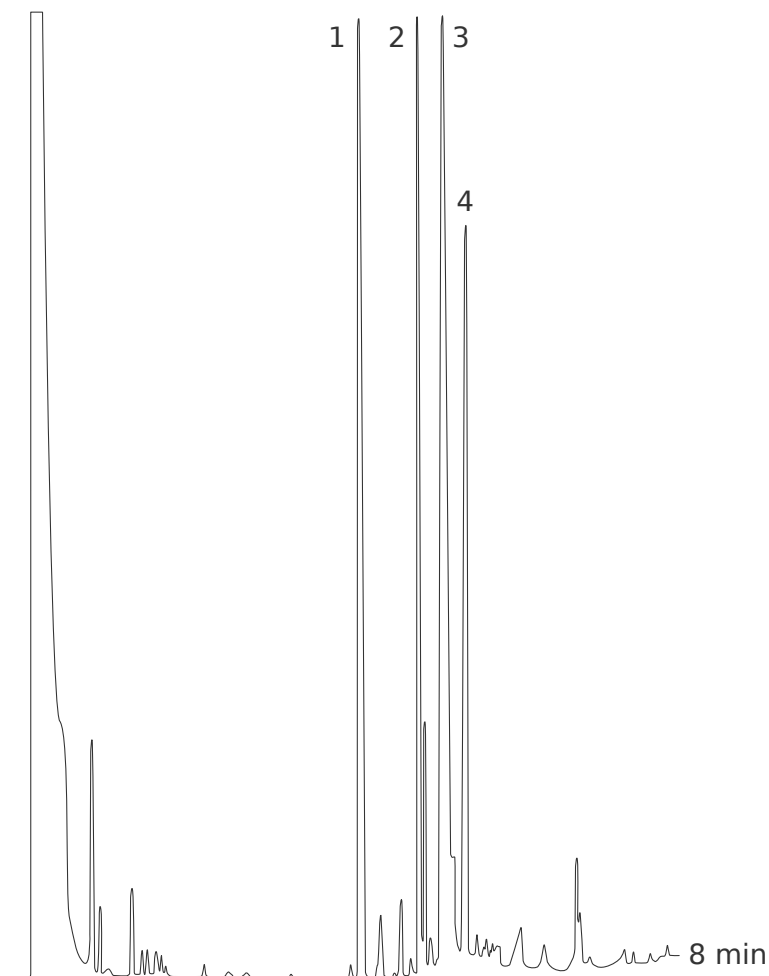




Drugs of Abuse - Cannabinoids

Column: **MEGA-1 - 0.25mm, 0.25 μ m, 15m**

Catalog Code: S-1-025-025-15

Conditions:

Injection: Split 280°C, 1.0 μ L, 1:50 Split Ratio.

Detector: FID, 320°C.

Oven Temp.: 120°C, 25°C/min, 310°C.

Carrier Gas: Hydrogen, 30kPa.

Peaks:

1. Docosane (I.S.)
2. CBD Cannabiol
3. THC Tetrahydrocannabiol
4. CBN Cannabiol

