

Faustime Essential Oil

(*Monocitrus australatica*, *Fortunella sp.*, *Citrus aurantifolia*)

Chiral Separation on MEGA-DEX Columns

Column: **MEGA-DEX DET Beta - 0.25mm, 0.25 μ m, 25m**

Catalog Code: DEX-DET-B-025-025-25

Conditions:


Injection: Split 250°C, 1 μ L, 1:50 Split Ratio.

Sample dilution: 1:200 in n-Hexane.

Contact us for additional information.

Peaks:

1. 1R 5R (+), 1R 5R (-) **Beta-Pinene**
2. 1R 5R (+), 1S 5S (-) **Sabinene**
3. 4S (-), 4R (+) **Limonene**
4. 3R (-), 3S (+) **Linalool**
5. 4S (-), 4R (+) **Alpha-Terpineol**

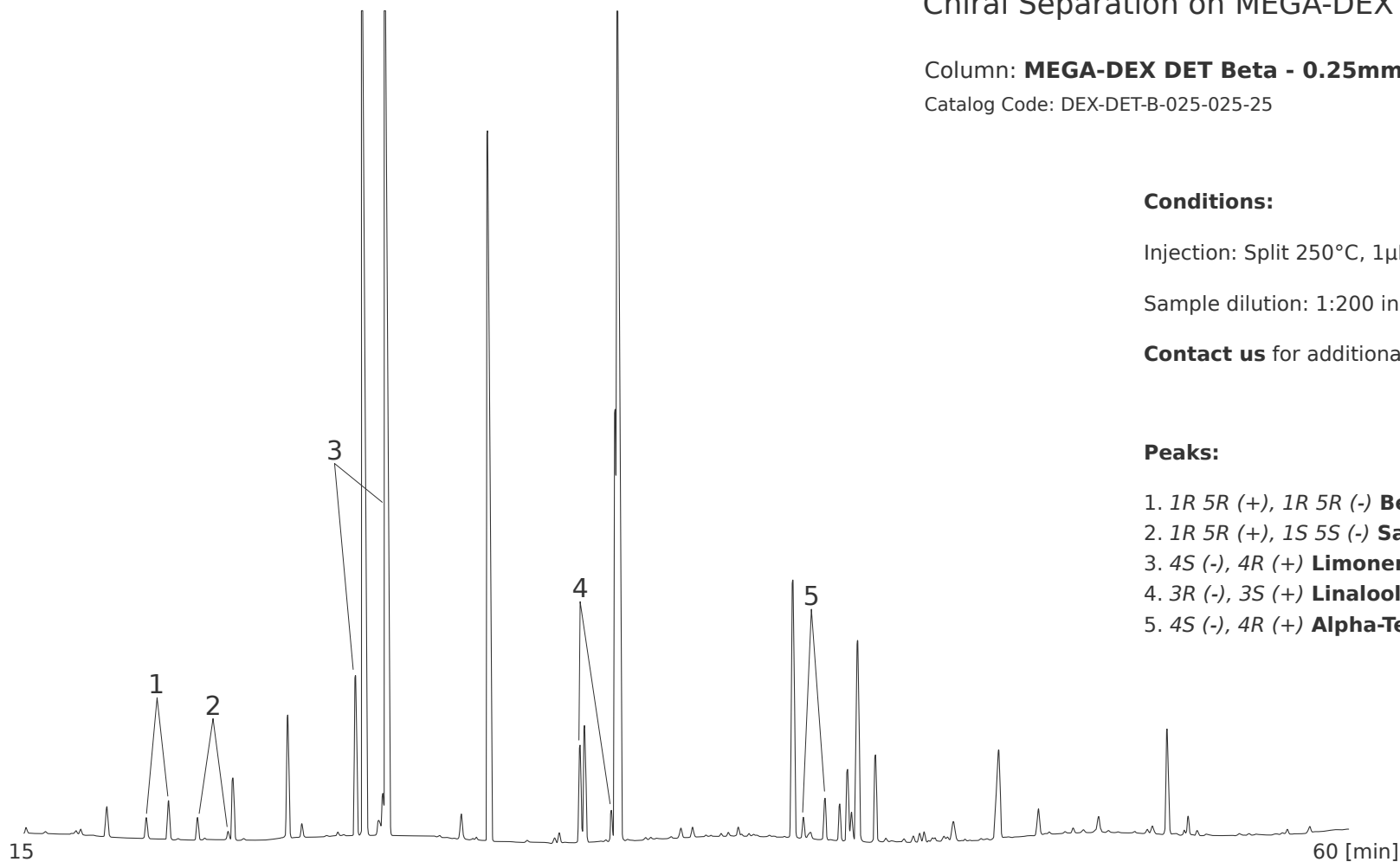

Acknowledgement: Prof. L.Mondello, P.Dugo et al. - Università di Messina, Dipartimento Farmaco-Chimico - Messina - Italy.

