


Hexobarbital (Pharmaceuticals) Chiral Separation on MEGA-DEX Columns

Column: **MEGA-DEX DMP Beta - 0.25mm, 0.25µm, 25m**

Catalog Code: DEX-DMP-B-025-025-25

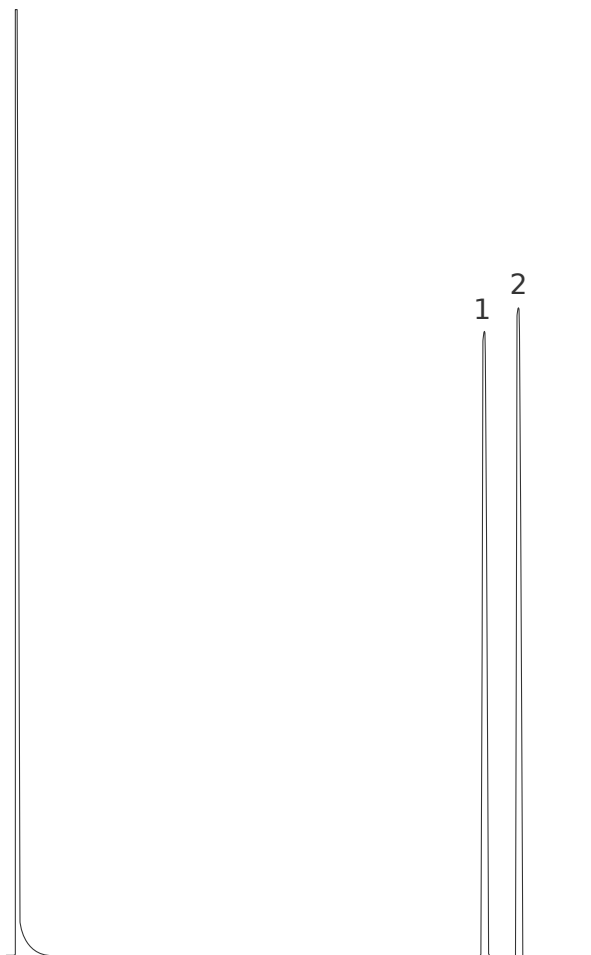

Conditions:

Injection: Split 250°C, 0.5µL, 1:50 Split Ratio.

Sample dilution: 1mg/1mL in n-Hexane.

Detector: FID 250°C.

Oven Temp.: 100°C, 2.5°C/min, 170°C

Carrier Gas: Hydrogen, 70kPa.

Peaks:

1. (+) Hexobarbital
2. (-) Hexobarbital

