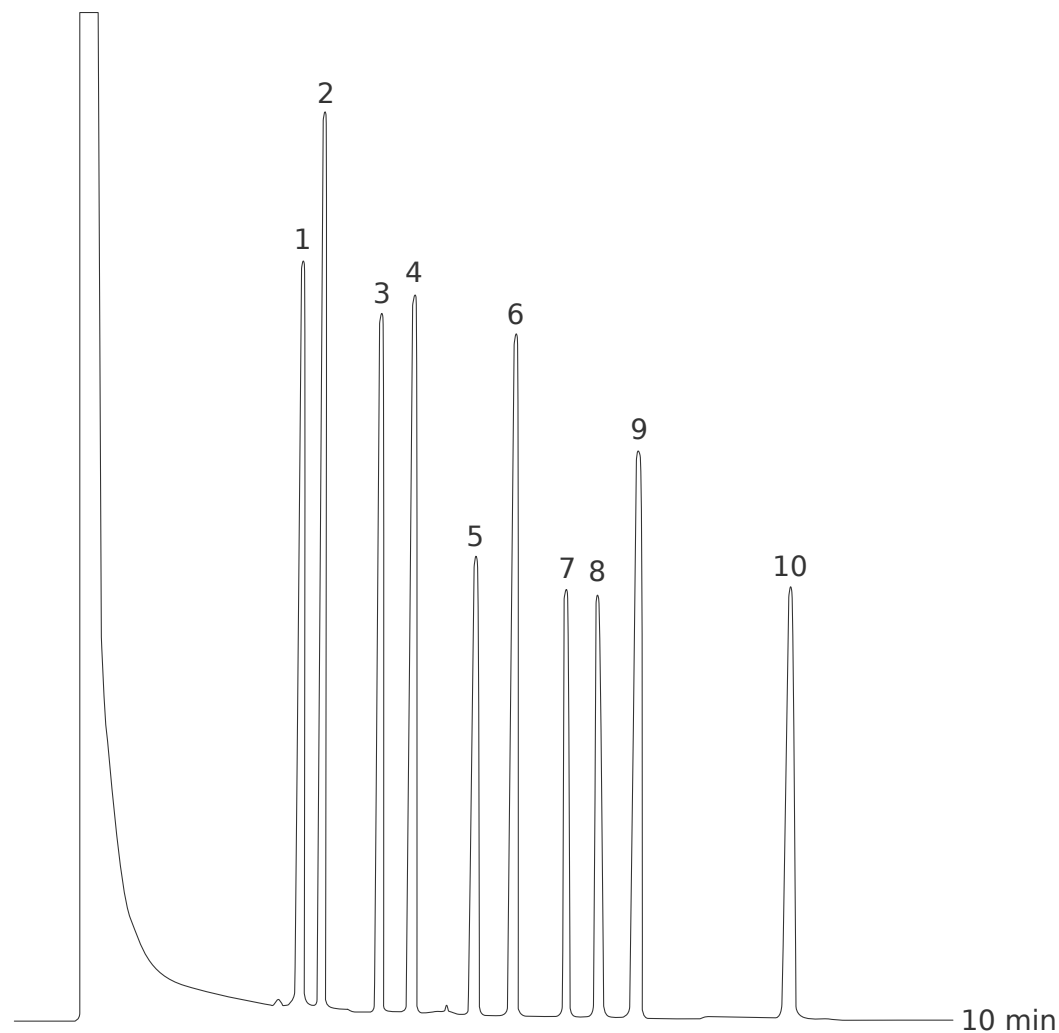




Nitrogen/Phosphorus Pesticides & Herbicides (Triazines) EPA Method 619

Column: **MEGA-WAX - 0.32mm, 0.25µm, 25m**

Catalog Code: S-WAX-032-025-25

Conditions:

Injection: Split 250°C, 1.5µL, 1:70 Split Ratio.

Detector: FID 300°C.

Oven Temp.: 230°C Isothermal.

Carrier Gas: Helium, 100 kPa.

Peaks:

CAS #:

1. Trietazine	1912-26-1
2. Prometon	1610-18-0
3. Propazine	139-40-2
4. Terbutylazine	5915-41-3
5. Atrazine	1912-24-9
6. Prometryne	7287-19-6
7. Terbutryne	886-50-0
8. Simazine	122-34-9
9. Ametryne	834-12-8
10. Simetryne	1014-70-6

