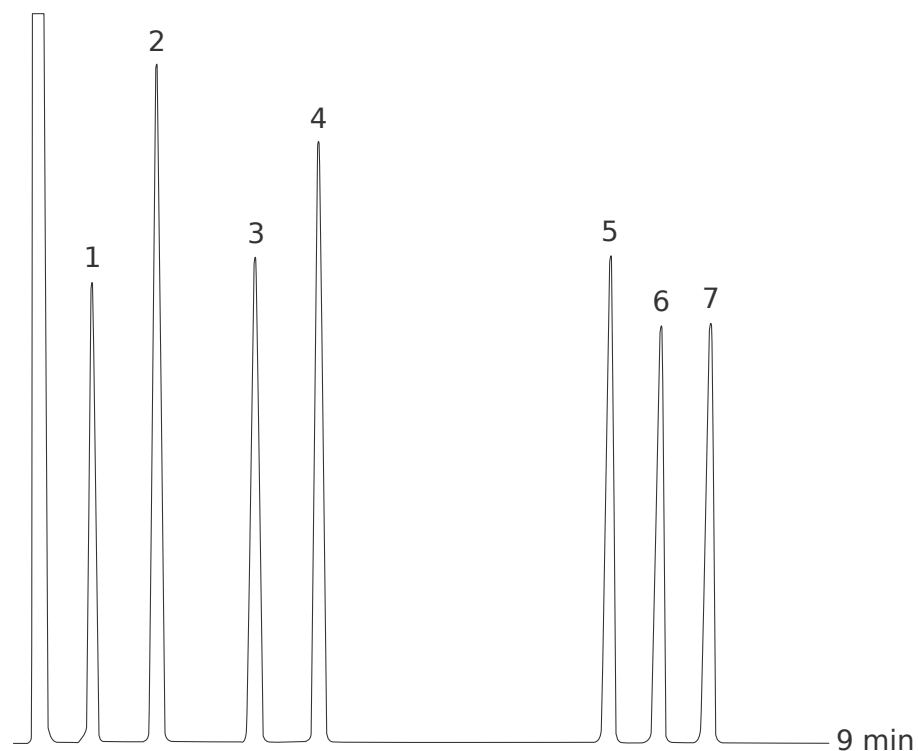




Purgeable Aromatics EPA Method 602

Column: **MEGA-WAX - 0.53mm, 1.00µm, 25m**

Catalog Code: S-WAX-053-100-25

Conditions:

Injection: Split 250°C, 1.5µL.

Sample: 2000µg/mL each in P&T Methanol.

Detector: FID 250°C.

Oven Temp.: 50°C, 8°C/min, 120°C.

Carrier Gas: Helium, 20 kPa.

Peaks:

1. Benzene
2. Toluene
3. Ethylbenzene
4. Chlorobenzene
5. 1,3-dichlorobenzene
6. 1,4-dichlorobenzene
7. 1,2-dichlorobenzene

