


Vodka on MEGA-WAX Spirit

Column: **MEGA-WAX Spirit, 0.32mm, 0.50µm, 50m**

Catalog Code: C-WAXSP-032-050-50

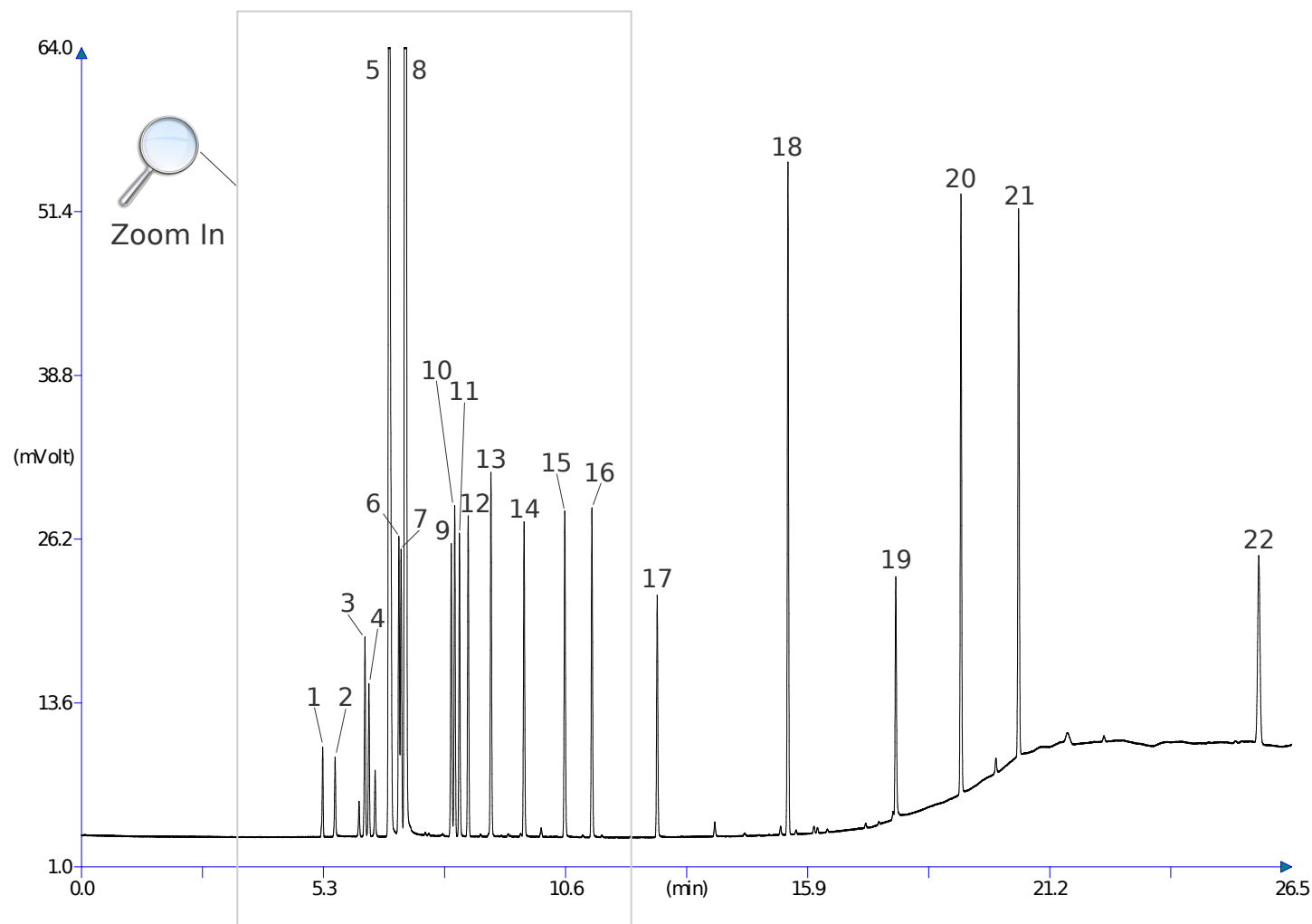

Conditions:

Injection: Split, 250°C, 1.0µL,
50mL/min Split Ratio.

Detector: FID, 250°C.

Oven Temp.: 60°C (2min), 10°C/min,
250°C.

Carrier Gas: Nitrogen, 0.8 mL/min.

Peaks:

1. Diethyl ether
2. Acetaldehyde
3. Acetone
4. Methylacetate
5. Methanol
6. MEK
7. 2-propanol
8. Ethanol
9. Isobutylacetate
10. 2-butanol
11. 1-propanol
12. Ethylbutirate
13. Crotonic aldehyde
14. Isobutanol
15. 1-butanol
16. Isopentanol
17. 1-pentanol
18. Hexanol
19. Benzaldehyde
20. Benzyl alcohol
21. Phenyl alcohol
22. Diethyl phthalate



